
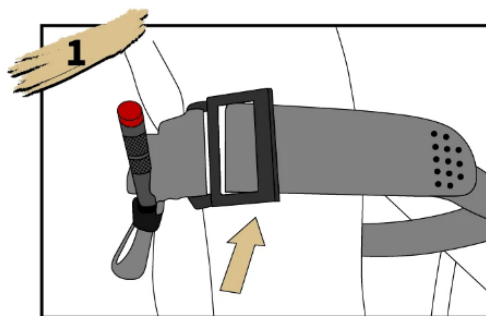
APLIKACE JEDNOU RUKOU

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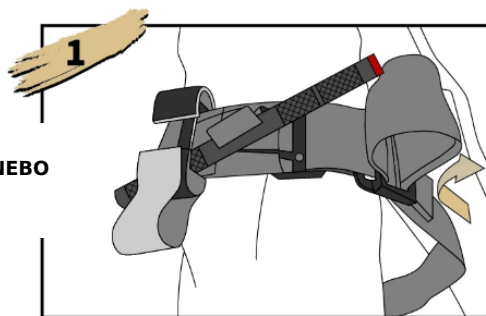
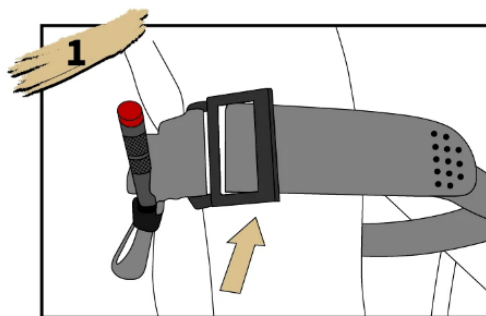
ONE-HANDED APPLICATION

Insert the injured limb through the loop in the strap and position the tourniquet 2 to 3 inches above the bleeding site, directly on the skin.

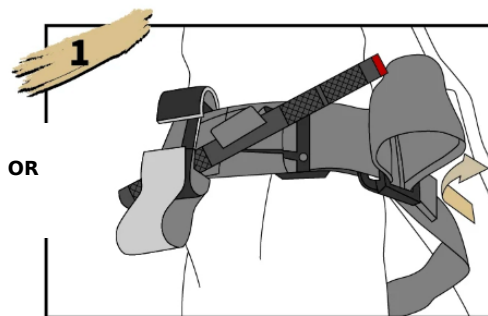OR


TWO-HANDED APPLICATION

Wrap the strap around the limb and pass the end through the openings in the buckle. Position the tourniquet 2 to 3 inches above the wound.


TOURNIQUET


INSTRUCTIONS FOR USE

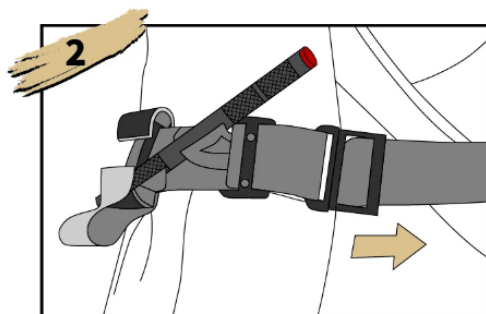
The strap must always pass around both the roller and the sliding teeth in the buckle. Pull the strap tight and secure it back onto itself around the limb, tight enough that three fingers cannot slide underneath. If possible, retighten and secure it.

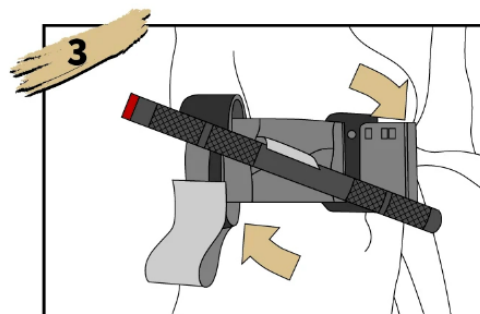Twist the windlass until the bleeding stops. Also check the distal pulse below the tourniquet near the wound; it should no longer be detectable.

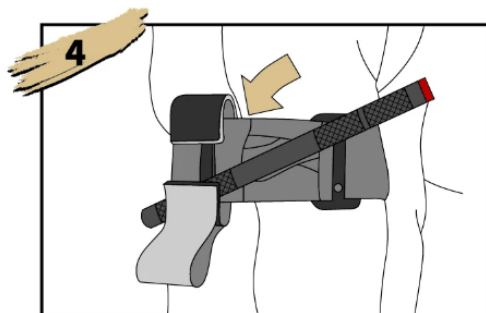Secure the windlass inside the clips. Check for bleeding and a distal pulse. If bleeding continues or a distal pulse is still present below the tourniquet, tighten it further or apply a second tourniquet above the first.

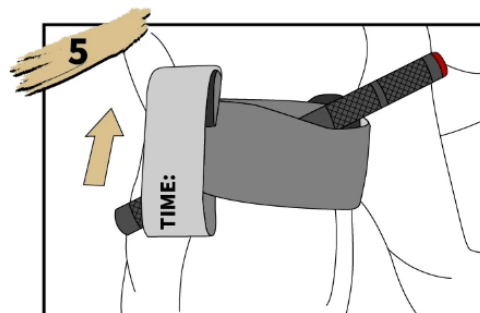Route the strap between the clips and over the windlass. Secure the windlass and strap with the TIME band. Record the application time in the marked area at the top. Consider additional securing with medical tape. Make sure the tourniquet remains visible to medical personnel.